

INCU-SHAKER™ 20L/LR

STACKABLE SHAKING INCUBATORS



Stackable design for high throughput culturing



Stainless steel chamber delivers precise temperature uniformity



MAGic-Clamp™ accessories for instant flask exchange



INCU-SHAKER™ 20L/LR STACKABLE SHAKING INCUBATOR



Convenient stacking

The included stacking brackets provide secure and easy mounting for stacking the Incu-Shaker 20L or 20LR, up to two units high. When stacked, the Incu-Shaker 20L and 20LR are still rated for continuous use with full flask capacity.

Stainless steel chamber

New to the Incu-Shaker series is the addition of a high grade stainless-steel chamber. The chamber is precision-engineered with consistent, repeatable culture growth in mind, and provides excellent temperature uniformity throughout the chamber.




Adjustable shelf

An included pull-out, stainless steel shelf allows for static incubation in addition to shaking incubation. Perfect for Petri dishes and other cell culture plates, the shelf can be adjusted for height as needed for additional capacity or flask height clearance.

Touch-screen controls

Shaking speed, time, and temperature can all be displayed and controlled on the full-color touch screen. The software also allows for more advanced programming options, including calibration, alarm deviations and the linking of consecutive programs.

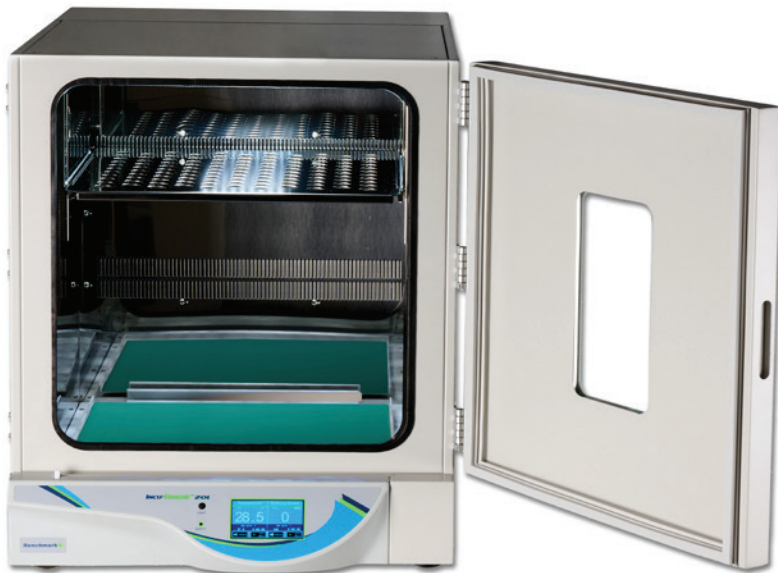

Technical Data:

Temperature Range (20L):	Ambient +5 to 60°C
Temperature Range (20LR):	4°C to 60°C
Temperature Accuracy:	±0.1°C at 37°C
Temperature Uniformity:	±0.5°C at 37°C
Temperature Increment:	0.1°C
Speed Range:	30 to 300 rpm
Speed Increment:	1 rpm
Mixing Orbit:	19mm (3/4 in.)
Time Range:	1 min. to 48 hours (or continuous)
Platform Dimensions:	18 x 18 in. (45 x 45 cm)
Chamber Dimensions: (W x D x H)	20 x 19.9 x 22.7 in. (51 x 50.5 x 57.7 cm)
Chamber Volume:	149L
Chamber Material:	Stainless Steel, 304
Exterior Dimensions: (W x D x H)	23.6 x 34.3 x 28.7 in. (60 x 87 x 73 cm)
Weight:	304 lbs / 138 kg
Electrical:	120V or 230V, 50-60Hz

Maximum Capacity Chart:

125ml Flasks	250ml Flasks	500ml Flasks	1000ml Flasks	2000ml Flasks	2.8L Fernbach	3000ml Flasks	5000ml Flasks
36	25	16	9	5	5	4	4



MAGic Clamp Platform
H1010-MR with Adjustable
MAGic Clamp (H1000-MR-CMB)

The Benchmark Incu-Shaker™ 20L and 20LR Stackable Shaking Incubators offer the highest capacity in the Incu-Shaker line. With a large shaking platform and a +22 inch chamber height, the 20L series is capable of high-throughput culture growth in flask sizes up to 5/6 liters. These shakers can be stacked (up to two units high), increasing the capacity while maximizing floor-space in the laboratory.


The internal chambers are constructed of high grade stainless-steel (304) and when combined with our proprietary heat distribution technology, the temperature is precisely maintained throughout the chamber. The powerful brushless induction motor drive is designed for years of trouble-free service, even under continuous use conditions with heavy loads. Noise and vibration are minimized by an advanced, computer designed motor isolation system that results in minimal vibration at all speed levels.

A convenient access port is included on the left side of both units, allowing for external temperature monitoring, or for power cables to run into the chamber without compromising the integrity of the environment. Additionally, each shaker includes an adjustable shelf that allows for static incubations to be performed simultaneously above the shaking platform.

A standard, rubber mat platform is included and is ideal for use with tissue culture flasks, Petri dishes, staining trays and other flat vessels. The convenient MAGic Clamp™ platform is also available for use with Erlenmeyer flasks and test tube racks. The unique, magnetic attachment method (PATENTED) is the easiest way to instantly change between flask clamps of different sizes.

Ordering Information:

H3020*	Incu-Shaker 20L with 18 x 18" non-slip rubber mat platform, 115V
H3022*	Incu-Shaker 20LR with 18 x 18" non-slip rubber mat platform, 115V
	<u>MAGic Clamp Accessories</u>
H1010-MR	MAGic Clamp Universal Platform, 18 x 18 inch
H1000-MR-CMB	MAGic Clamp Adjustable Flask Clamp, for 125, 250 & 500ml Erlenmeyer flasks
H1000-MR-50	Magic Clamp Flask Clamp, for 50ml Erlenmeyer flask (max. 49)
H1000-MR-125	Magic Clamp Flask Clamp, for 125ml Erlenmeyer flask (max. 36)
H1000-MR-250	Magic Clamp Flask Clamp, for 250ml Erlenmeyer flask (max. 25)
H1000-MR-500	Magic Clamp Flask Clamp, for 500ml Erlenmeyer flask (max. 16)
H1000-MR-1000	Magic Clamp Flask Clamp, for 1L Erlenmeyer flask (max. 9)
H1000-MR-2000	Magic Clamp Flask Clamp, for 2L Erlenmeyer flask (max. 5)
H1000-MR-3000	Magic Clamp Flask Clamp, for 3L Erlenmeyer flask or 2.8L Fernbach (max. 4)
H1000-MR-5000	Magic Clamp Flask Clamp, for 5L Erlenmeyer/Fernbach flask (max. 4)
H1000-MR-MP	Magic Clamp Microplate Clamp, for 1 microplate or deep-well plate
H1000-MR-T15	Magic Clamp Tube Rack, 32 x 15ml
H1000-MR-T50	Magic Clamp Tube Rack, 15 x 50 ml
H1000-MR-1550	Magic Clamp Tube Rack combination rack, 20 x 50 ml and 30 x 15 ml
H1000-MR-T600	Magic Clamp Tube Rack, 4 x 500/600ml bottles
H1000-MR-TSR	Magic Clamp Slanted Tube Rack Holder, adjustable angle
	*115V with USA Plug. For 230V with EU plug, add (-E)

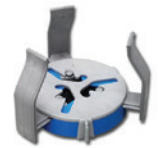

H1000-MR-CMB




H1000-MR-1550



H1000-MR-MP



H1000-MR-T600

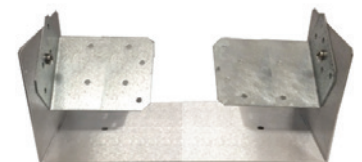

H1000-MR-TSR

Ordering Information (Flasks):

F4060-F	125 ml Erlenmeyer flask, PETG, flat base, disposable, qty 24
F4061-F	250 ml Erlenmeyer flask, PETG, flat base, disposable, qty 12
F4062-F	500 ml Erlenmeyer flask, PETG, flat base, disposable, qty 12
F4063-F	1000 ml Erlenmeyer flask, PETG, flat base, disposable, qty 6
F4060-B	125 ml Erlenmeyer flask, PC, baffled base, autoclavable/reusable, qty 24
F4061-B	250 ml Erlenmeyer flask, PC, baffled base, autoclavable/reusable, qty 12
F4062-B	500 ml Erlenmeyer flask, PC, baffled base, autoclavable/reusable, qty 12
F4063-B	1000 ml Erlenmeyer flask, PC, baffled base, autoclavable/reusable, qty 6
F4086	2.8L Jumbo Fernbach flask, styrene/butadiene, baffled base, qty 6
F4087	5L Jumbo Fernbach flask, styrene/butadiene, baffled base, qty 4



Also Available:



INCU-SHAKER™
Mini

Compact footprint with exceptional performance and reliability



INCU-SHAKER^{CO2}
Mini

Compact CO₂ incubator with integrated shaker